


Disposable Micro Vascular Clip

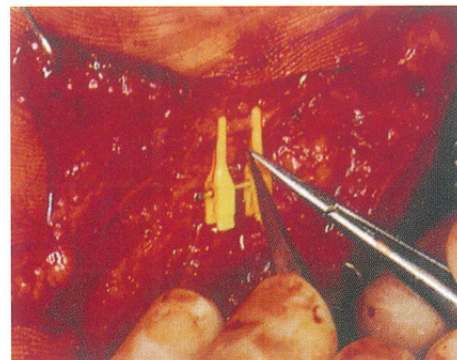

TK Series

For Artery : Yellow


Single type : 4 sizes

Double type : 4 sizes

For Vein : Green

Single type : 4 sizes


Double type : 4 sizes

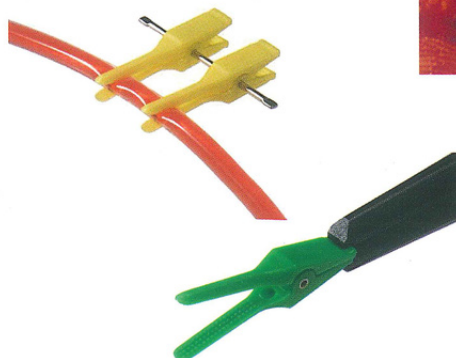

※used in case of several operations especially micro surgery by prastic surgeons.
(Flap transplatation, Blood vessel or lymph anastomosis, Forming dialysis shunt, etc...)


Green : for vein Yellow : for artery		Product code		Clipping power	Clip length	Proper vessel diameter	Quantity/box
		Single clip	Double clip				
		TKF-1-15	TKF-2-15	0.15N(15gf)	11mm	0.2 ~ 1mm	10pcs
		TKF-1-30	TKF-2-30	0.29N(30gf)			
		TKS-1-20	TKS-2-20	0.20N(20gf)	14.5mm	MAX 1mm	10pcs
		TKS-1-40	TKS-2-40	0.39N(40gf)			
		TKM-1-30	TKM-2-30	0.29N(30gf)	17mm	1 ~ 2mm	10pcs
		TKM-1-60	TKM-2-60	0.59N(60gf)			
		TKL-1-60	TKL-2-60	0.59N(60gf)	21.5mm	2 ~ 4mm	10pcs
		TKL-1-120	TKL-2-120	1.18N(120gf)			

Holding Mechanism
(Frontal cut section)

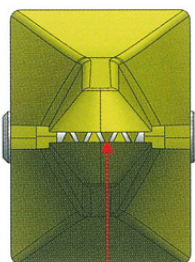

Arrangement of fine theeth on the upper and lower jaw surgaces can bitten alternately

Holding Mechanism
(Lateral view)

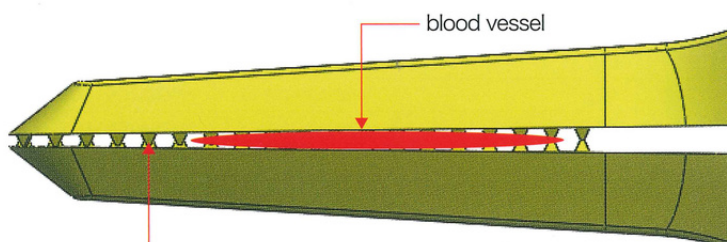

Pallarel holding power on the blood vessels can be testified with suitable diameter of the vessels